


Medical Oxygen Regulator



Nova 40

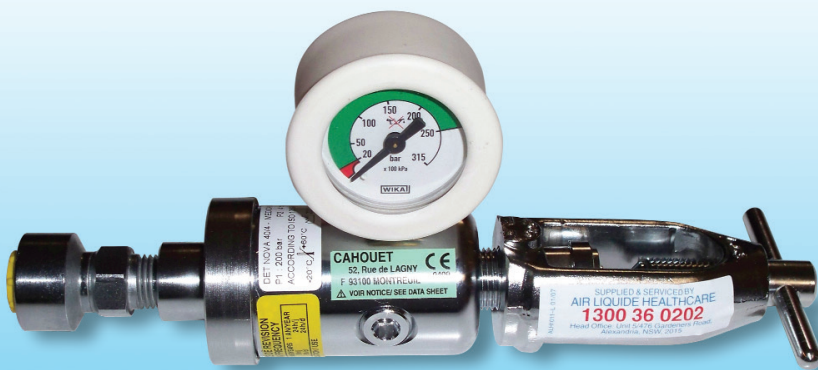

- Simple to use
- Low Maintenance
- Maximum Safety

With its simple design and superior engineering, the Nova-40 regulator is a safer and more cost effective¹ solution for medical oxygen cylinders.


Key features include:

- Full brass construction for complete safety at high pressure
- A built-in safety relief valve
- 5-yearly service interval²
- 360° rotatable body for optimum positioning
- Compliant with AS3840.1 (1998)
- Included on the ARTG as a medical device

¹ In comparison to ordinary regulators

² Depending on frequency of use

NOVA-40 Regulator

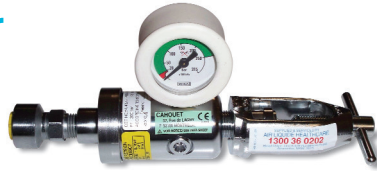

NOVA-40 regulator is applicable wherever the administration of medical oxygen to a patient is required from a pressurised gas cylinder. It is ideal for use in the delivery of oxygen in emergency situations, dental and home oxygen therapy applications.

Specification

- Rated inlet pressure: 20,000kPa
- Outlet pressure: 400kPa
- Flow rate: 50 Litre/min
- Operating temperature: -20°C to +60°C
- Weight: 960 grams
- Body: Chrome plated brass

Standard Configuration

- Pin index yoke for cylinder connection
- Sleeve index outlet (S.I.S.)
- Color coded content gauge with a rubber guard
- Sintered bronze inlet filter
- Safety relief valve
- Each regulator is serialized

Note: S.I.S as per AS2896 (1998)

Low Maintenance Regulator = Operational Savings and Peace of Mind

Whilst ordinary regulators need a 3 year overhaul service and yearly filter change, the Nova 40 only requires a 5 year overhaul service and no filter change under normal usage*. This means that you save in excess of \$200 in the first 5 years!

** Savings may vary due to different types of regulators used in the comparison; the savings estimated above is based on normal usage of 8 hours a day.*

AIR LIQUIDE Healthcare Equipment Maintenance Program

We offer a Preventative Maintenance program for all Medical Gas Equipment including regulators, wall outlets and gas manifolds. Please contact us for more information.

Contact us now to know more about our services on 1300 36 02 02.

NSW State Office

Unit 5/476 Gardeners Road
ALEXANDRIA NSW 2015
Tel: 02 9364 7474
Fax: 02 8338 9797

QLD State Office

25 Campbell Street
BOWEN HILLS QLD 4006
Tel: 07 3252 8039
Fax: 07 3252 4850

VIC & TAS State Office

40 Bunnet Street
NORTH SUNSHINE VIC 3020
Tel: 03 9310 1200
Fax: 02 9364 7482

SA State Office

95 Fullarton Road
KENT TOWN SA 5067
Tel: 08 8331 1648
Fax: 08 8331 1849

AIR LIQUIDE Healthcare Pty. Ltd.
<http://www.airliquidehealthcare.com.au>
Toll free 1300 36 02 02

