

Matrix Scavenger



Part # 3570S

**Comfortable
patient circuit
with effective
scavenging
design**

Matrix Scavenger is a fully autoclavable patient circuit for nitrous oxide sedation suitable for use in dental surgeries.

- Efficient gas delivery and scavenging
- Fully autoclavable
- Comfortable for the patient

BENEFITS OF MATRX SCAVENGER

✓ Better performance

- Lightweight design provides better balance on the patient's face. The result is a better seal around the nose for optimal gas delivery and scavenging.
- Unique exterior scavenging cone design allows removal of peripheral waste gases.

✓ Safe to use

- Fully autoclavable
- Latex-free

✓ Fewer parts... easier to assemble

- Simple 8-piece circuit for easy assembly

✓ More comfortable for your patient

- Weighs 50% less than other units on the market - your patients are more comfortable and relaxed.

Did you know?

Scavenging waste gases during the administration of nitrous oxide sedation is an important requirement listed in the "Guidelines on Conscious Sedation for Dental Procedures" (PS21) published by ANZCA*.

* Australian and New Zealand College of Anaesthetists



Matrix Scavenger only 8 parts to assemble

NASAL HOODS

See (1) figure 1. below

• Reusable nasal hoods

- Fully autoclavable
- Comes in 3 sizes



Small: Part # 848S
Medium: Part # 847S
Large: Part # 846S

• Dynamite™

Disposable, scented nasal hoods

- 5 delicious scents
- Amusing dinosaur stickers
- Latex-free



Small: Part # 1546S
Medium: Part # 1547S
Large: Part # 1548S

REPLACEMENT PARTS

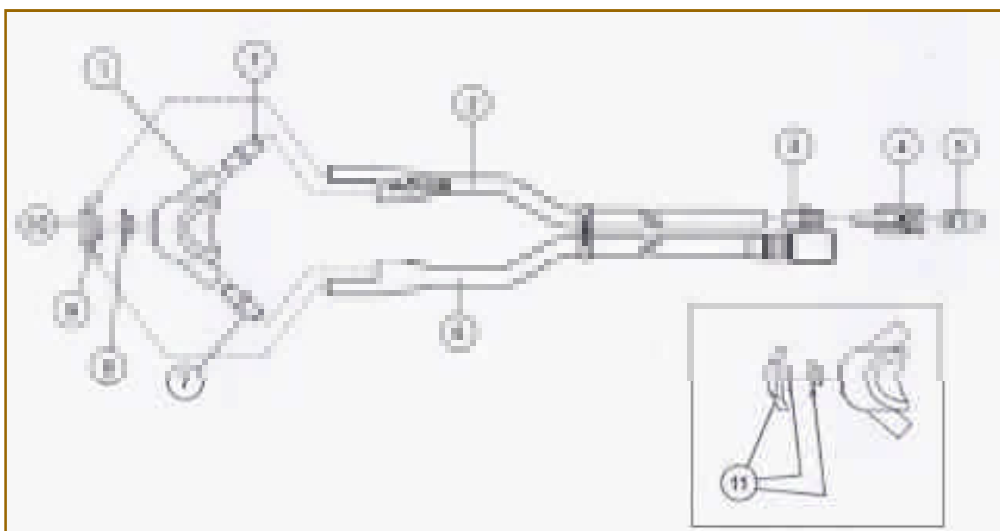


Figure 1. Matrix Scavenger parts

Key (figure 1)	Replacement parts	Part Number
2	Vacuum Tubing	3772S
3	Vacuum Nipple (small orifice)	1279S
4	Shut-Off Valve	1272S
5	Tapered Vacuum Adaptor *	1448S
6	Gas Tubing w/slide	3771S
7	Hood/Hose Connector	1075S
8	Cage Valve	3773S
9	Non-Breathing Valve (5pk)	3774S
10	Cone Assembly	3776S
11	Cone only	3775S



For further information, please contact your AIR LIQUIDE Healthcare representative:
Toll free 1300 36 02 02 / <http://www.airliquidehealthcare.com.au>



Quality Endorsed Company
ISO 9001 Lic: 1457601
Standards Australia