


INVACARE®  
**PERFECTO<sub>2</sub>™ V**  
OXYGEN CONCENTRATOR

---

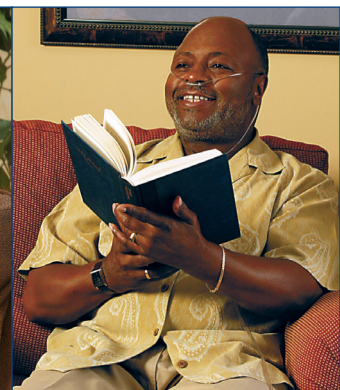

---

The **Invacare® Perfecto<sub>2</sub>™ V Oxygen Concentrator** is an economical addition to the premier Invacare stationary concentrator line. This 5-liter concentrator is small, quiet and lightweight and comes standard with the **Invacare® SensO<sub>2</sub>® Oxygen Sensor**. The Perfecto<sub>2</sub> V unit incorporates many of the best features of the Invacare® Perfecto<sub>2</sub>™ Concentrator, including its contemporary, patient preferred external design and color.

---

## Features

Pressure-based system is superior to timing-based systems  
 Self diagnostic electronics minimize troubleshooting time  
 Top handle for ease of mobility  
 Battery-free power loss alarm, compressor 35 psi pressure-relief valve and low flow alarm

## Provider Benefits

Robust design comes standard with Invacare® SensO2® oxygen monitor to simplify preventative maintenance schedules and lower costs  
 All Perfecto2 V stationary concentrators include a HomeFill system port for inventory flexibility  
 Patient-accessible reset button helps minimize emergency service calls

## Patient Benefits

User friendly design with easy to see control panel, humidifier bottle holder and filter access  
 Low noise level and neutral color help to blend into any home environment  
 Reliable concentrator provides unlimited, worry-free refillable oxygen supply

## HCPCS Code: E1390



1. Invacare® HomeFill® Oxygen System connection
2. HomeFill compressor and Ready Rack have been updated to coordinate with the Perfecto<sub>2</sub> Concentrator. Use with Invacare cylinders
3. Easy to maintain filter (not pictured)



**Invacare Corporation**  
 www.invacare.com

### USA

One Invacare Way  
 Elyria, Ohio  
 44035-4190  
 (800) 333-6900

### Canada

570 Matheson Blvd. E.,  
 Unit 8  
 Mississauga, Ontario  
 L4Z 4G4 Canada  
 (800) 668-5324

©2017

Invacare Corporation.  
 All rights reserved.  
 Trademarks are identified by the symbols™ and ©. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Specifications are subject to change without notification.  
 Form No. 09-064  
 Rev. 06/17  
 170583

Model #	Dimensions	Weight	Shipping Weight	Quantity	Oxygen Flow	Flow Range	Maximum Outlet Pressure	Pressure Relief Mechanism
IRC5PO2V	23" H 13" W 11.5" D	39 lb.	45 lb.	1 Per Case 27 Per Pallet	5 LPM	0.5-5 LPM	5+/- 0.5 PSI	35 PSI +/- 5 PSI
	Oxygen Purity	Filter	Electrical Requirements	Rated Current Input	Average Power Consumption	Sound Level	Altitude	Limited Warranty
	95.6%-87% At All Flow Rates	Cabinet Compressor Inlet Outlet HEPA	120 VAC +/- 10%, 60 Hz	3.0 A	320 W @ 5 LPM	48 db(A)	8,000 ft.	3 Years

All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further detail.

