

Safety Data Sheet

MEDICAL CARBON DIOXIDE, Compressed & Liquefied Gas

Date of first issue: 27/09/2010 SDS reference: ALH062 Revised date: 13/11/2017

Supersedes: 16/11/2012

H280

Version: 5.0

Warning



SECTION 1: Identification of the substance/mixture and of the company/undertaking

1.1. Product identifier	
Trade name	: Medical Carbon dioxidem EP Grade
SDS no	: ALH062
Chemical description	: Carbon dioxide
	CAS-No. : 124-38-9
Chemical formula	: CO2
1.2. Relevant identified uses of the	substance or mixture and uses advised against
Relevant identified uses	 Medical: Used with medical oxygen to stimulate breathing response. Insufflation has for minimal invasive surgery. Cryotherapy. Cryopreservation.
	Pulp sensitivity testing.
Uses advised against	: Consumer use.
1.3. Details of the supplier of the sa	fety data sheet
Company identification	 Air Liquide Healthcare Pty Limited Unit 5, 476 Gardeners Road 2015 Alexandria NSW Australia +61 2 9364 7474 ALHEnguiries@AirLiguide.com
1.4. Emergency telephone number	
Emergency telephone number	: 1800 812 588

SECTION 2: Hazards identification

2.1. Classification of the substance or mixture

Classification according to WHS Regulation

Physical hazards Gases under pressure : Liquefied gas

2.2. Label elements

Classification according to WHS Regulation

Hazard pictograms

	GHS04	
Signal word	: Warning	
Hazard statements	: H280 - Contains gas under p	ressure; may explode if heated
Air Liquide Healthcare Pty Limited	EN (English)	SDS Ref.: ALH062



SDS Ref.: ALH062

Precautionary statements

- Storage : P403 - Store in a well-ventilated place..

2.3. Other hazards

Asphyxiant in high concentrations.
 Contact with liquid may cause cold burns/frostbite.
 In high concentrations CO2 cause rapid circulatory insufficiency. Symptoms are headache, nausea and vomiting, which may lead to unconsciousness.

SECTION 3: Composition/information on ingredients

3.1. Substances

Name	Product identifier	%	Classification according to WHS Regulation
Carbon dioxide	(CAS-No.) 124-38-9	100	Press. Gas (Liq.), H280

Contains no other components or impurities which will influence the classification of the product.

Full text of R-phrases see section 16. Full text of H-statements see section 16.

3.2. Mixtures : Not applicable

SECTION 4: First aid measures

Description of first aid measures 4.1. - Inhalation : Remove victim to uncontaminated area wearing self contained breathing apparatus. Keep victim warm and rested. Call a doctor. Perform cardiopulmonary resuscitation if breathing stopped. - Skin contact : In case of frostbite spray with water for at least 15 minutes. Apply a sterile dressing. Obtain medical assistance. - Eye contact : Immediately flush eyes thoroughly with water for at least 15 minutes. - Ingestion : Ingestion is not considered a potential route of exposure. <u>4.2.</u> Most important symptoms and effects, both acute and delayed : In high concentrations may cause asphyxiation. Symptoms may include loss of mobility/consciousness. Victim may not be aware of asphyxiation. Low concentrations of CO2 cause increased respiration and headache. Refer to section 11.

4.3. Indication of any immediate medical attention and special treatment needed

: None.

SECTION 5: Firefighting measures

- Suitable extinguishing media	: Water spray or fog.
- Unsuitable extinguishing media	Do not use water jet to extinguish.
5.2. Special hazards arising from	the substance or mixture
Specific hazards	: Exposure to fire may cause containers to rupture/explode.
Hazardous combustion products	: None.
5.3. Advice for fire-fighters	



SDS Ref.: ALH062

Specific methods	 Use fire control measures appropriate for the surrounding fire. Exposure to fire and heat radiation may cause gas receptacles to rupture. Cool endangered receptacles with water spray jet from a protected position. Prevent water used in emergency cases from entering sewers and drainage systems. If possible, stop flow of product. Use water spray or fog to knock down fire fumes if possible. Move containers away from the fire area if this can be done without risk.
Special protective equipment for fire fighters	 In confined space use self-contained breathing apparatus. Standard protective clothing and equipment (Self Contained Breathing Apparatus) for fire fighters. Standard EN 137 - Self-contained open-circuit compressed air breathing apparatus with full face mask. Standard EN 469 - Protective clothing for firefighters. Standard - EN 659: Protective gloves for firefighters.
Hazchemcode	: 2T

SECTION 6: Accidental release measures

6.1. Personal precautions, protective equipment and emergency procedures

		: Try to stop release.
		Evacuate area.
		Wear self-contained breathing apparatus when entering area unless atmosphere is proved to be safe.
		Ensure adequate air ventilation.
		Prevent from entering sewers, basements and workpits, or any place where its accumulation can be dangerous.
		Act in accordance with local emergency plan.
		Stay upwind.
		Oxygen detectors should be used when asphyxiating gases may be released.
<u>6.2.</u>	Environmental precautions	
		: Try to stop release.
<u>6.3.</u>	Methods and material for containme	nt and cleaning up
		: Keep area evacuated and free from ignition sources until any spilled liquid has evaporated (ground free from frost).
6.4.	Reference to other sections	
		: See also sections 8 and 13.

SECTION 7: Handling and storage

7.1. Precautions for safe handling

Do not smoke while handling product. Use only properly specified equipment which is suitable for this product, its supply pressure and temperature. Contact your gas supplier if in doubt. Avoid suck back of water, acid and alkalis. Do not breathe gas. Avoid release of product into atmosphere. Containers, which contain or have contained flammable or explosive substances, must not be inerted with liquid carbon dioxide. Potential production of solid CO2 particles must be ruled out. In order to rule out potential electrostatic discharge production, the system must be adequately grounded.		Use only properly specified equipment which is suitable for this product, its supply pressure and temperature. Contact your gas supplier if in doubt. Avoid suck back of water, acid and alkalis. Do not breathe gas. Avoid release of product into atmosphere. Containers, which contain or have contained flammable or explosive substances, must not be
--	--	--

Air Liquide H E A L T H C A R E	MEDICAL CARBON DIOXIDE, Compressed & Liquefied Gas
	SDS Ref.: ALH062
Safe handling of the gas receptacle	 Refer to supplier's container handling instructions. Do not allow backfeed into the container. Protect cylinders from physical damage; do not drag, roll, slide or drop. When moving cylinders, even for short distances, use a cart (trolley, hand truck, etc.) designed to transport cylinders. Leave valve protection caps in place until the container has been secured against either a wall or bench or placed in a container stand and is ready for use. If user experiences any difficulty operating cylinder valve discontinue use and contact supplier. Never attempt to repair or modify container valves or safety relief devices. Damaged valves should be reported immediately to the supplier. Keep container valve outlets clean and free from contaminants particularly oil and water. Replace valve outlet caps or plugs and container caps where supplied as soon as container is disconnected from equipment. Never attempt to transfer gases from one cylinder/container to another. Never use direct flame or electrical heating devices to raise the pressure of a container. Do not remove or deface labels provided by the supplier for the identification of the cylinder contents. Suck back of water into the container must be prevented. Open valve slowly to avoid pressure shock.
7.2. Conditions for safe storage, in	ncluding any incompatibilities
	 Observe all regulations and local requirements regarding storage of containers. Containers should not be stored in conditions likely to encourage corrosion. Container valve guards or caps should be in place. Containers should be stored in the vertical position and properly secured to prevent them from falling over. Stored containers should be periodically checked for general condition and leakage. Keep container below 50°C in a well ventilated place. Store containers in location free from fire risk and away from sources of heat and ignition. Keep away from combustible materials.
7.3. Specific end use(s)	

Specific end use(s)

: None.

SECTION 8: Exposure controls/personal protection

8.1. **Control parameters**

Carbon dioxide (124	-38-9)		
OEL : Occupational E	xposure Limits		
Australia	TWA (mg/m³)	9000 mg/m ³	
	TWA (ppm)	5000 ppm	
	STEL (mg/m ³)	54000 mg/m ³	
	STEL (ppm)	30000 ppm	

DNEL (Derived-No Effect Level) : No data available.

PNEC (Predicted No-Effect Concentration) : No data available.

8.2. Exposure controls

8.2.1. Appropriate engineering controls

Provide adequate general	and local exhaust ventilation.

Systems under pressure should be regularly checked for leakages. Ensure exposure is below occupational exposure limits (where available). Oxygen detectors should be used when asphyxiating gases may be released. Consider the use of a work permit system e.g. for maintenance activities. CO2 detectors should be used when CO2 may be released.

8.2.2. Individual protection measures, e.g. personal protective equipment

: A risk assessment should be conducted and documented in each work area to assess the risks related to the use of the product and to select the PPE that matches the relevant risk. The following recommendations should be considered:

PPE compliant to the recommended Australian or EN/ISO standards should be selected.



SDS Ref.: ALH062

Eye/face protection	: Wear goggles when breaking transfer connections. Standard EN 166 - Personal eye-protection – specifications, AS/NZS 1337.1 – Eye and face protectors for occupational applications
Skin protection	
- Hand protection	 Wear working gloves when handling gas containers. Standard EN 388 - Protective gloves against mechanical risk; AS/NZS – 2161.3 – Occupational protective gloves, protection against mechanical risks. Wear cold insulating gloves when breaking transfer connections. Standard EN 511 - Cold insulating gloves.
- Other	 Wear safety shoes while handling containers. Standard EN ISO 20345 - Personal protective equipment - Safety footwear; AS/NZS 2210.3 – Occupational protective footwear, specifications for safety footwear.
Respiratory protection	 Gas filters may be used if all surrounding conditions e.g. type and concentration of the contaminant(s) and duration of use are known. Use gas filters with full face mask, where exposure limits may be exceeded for a short-term period, e.g. connecting or disconnecting containers. Gas filters do not protect against oxygen deficiency. Self contained breathing apparatus (SCBA) or positive pressure airline with mask are to be used in oxygen-deficient atmospheres. Standard EN 14387 - Gas filter(s), combined filter(s) and full face mask. Standard EN 137 - Self-contained open-circuit compressed air breathing apparatus with full face mask.
Thermal hazards	: None in addition to the above sections
8.2.3. Environmental exposure controls	

: None necessary.

SECTION 9: Physical and chemical properties

9.1. Information on basic physical and chemical properties

Appearance

 Physical state at 20°C / 101.3kPa 	: Gas.
Colour	: Colourless.
Odour	: No odour warning properties.
Odour threshold	: Odour threshold is subjective and inadequate to warn of overexposure.
pH value	: Not applicable for gases and gas mixtures.
Molar mass	: 44 g/mol
Melting point	: -78.5 °C At atmospheric pressure dry ice sublimes into gaseous carbon dioxide.
Boiling point	: -56.6 °C
Flash point	: Not applicable for gases and gas mixtures.
Critical temperature [°C]	: 30 °C
Evaporation rate (ether=1)	: Not applicable for gases and gas mixtures.
Flammability range	: Non flammable.
Vapour pressure [20°C]	: 57.3 bar(a)
Vapour pressure [50°C]	: Not applicable.
Relative density, gas (air=1)	: 1.52
Relative density, liquid (water=1)	: 0.82
Solubility in water	: 2000 mg/l Completely soluble.



SDS Ref.: ALH062

Partition coefficient n-octanol/water [log Kow]	: 0.83
Auto-ignition temperature	: Non flammable.
Decomposition point [°C]	: Not applicable.
Viscosity [20°C]	: No reliable data available.
Explosive Properties	: Not applicable.
Oxidising Properties	: Not applicable.
9.2. Other information	
Other data	: Gas/vapour heavier than air. May accumulate in confined spaces, particularly at or below ground level.

SECT	ON 10: Stability and reactivity	
10.1.	Reactivity	
		: No reactivity hazard other than the effects described in sub-sections below.
<u>10.2.</u>	Chemical stability	: Stable under normal conditions.
<u>10.3.</u>	Possibility of hazardous reactions	
<u>10.4.</u>	Conditions to avoid	: None.
10.4.	Conditions to avoid	: Avoid moisture in installation systems.
<u>10.5.</u>	Incompatible materials	
		: None. For additional information on compatibility refer to ISO 11114.
<u>10.6.</u>	Hazardous decomposition products	: None.

SECTION 11: Toxicological information

11.1. Information on toxicological effects

Acute toxicity	: Unlike simple asphyxiants, carbon dioxide has the ability to cause death even when normal oxygen levels (20-21%) are maintained. 5% CO2 has been found to act synergistically to increase the toxicity of certain other gases (CO, NO2). CO2 has been shown to enhance the production of carboxy- or met-hemoglobin by these gases possibly due to carbon dioxide's stimulatory effects on the respiratory and circulatory systems. For more information, see 'EIGA Safety Info 24: Carbon Dioxide, Physiological Hazards' at www.eiga.eu.
Skin corrosion/irritation	: No known effects from this product.
Serious eye damage/irritation	: No known effects from this product.
Respiratory or skin sensitisation	: No known effects from this product.
Germ cell mutagenicity	: No known effects from this product.
Carcinogenicity	: No known effects from this product.
Toxic for reproduction : Fertility	: No known effects from this product.
Toxic for reproduction : unborn child	: No known effects from this product.
STOT-single exposure	: No known effects from this product.
STOT-repeated exposure	: No known effects from this product.
Aspiration hazard	: Not applicable for gases and gas mixtures.



SDS Ref.: ALH062

SECTION 12: Ecological information	on
12.1. Toxicity	
Assessment	: No ecological damage caused by this product.
12.2. Persistence and degradability	
Assessment	: No ecological damage caused by this product.
12.3. Bioaccumulative potential	
Assessment	: No ecological damage caused by this product.
12.4. Mobility in soil	
Assessment	: Because of its high volatility, the product is unlikely to cause ground or water pollution. Partition into soil is unlikely.
12.5. Results of PBT and vPvB assessr	nent
Assessment	: No data available. Not classified as PBT or vPvB.
12.6. Other adverse effects	
	: No known effects from this product.
Effect on the ozone layer	: None.
Global warming potential [CO2=1]	: 1
Effect on global warming	: Contains greenhouse gas(es). When discharged in large quantities may contribute to the greenhouse effect.
SECTION 13: Disposal considerati	ons
12.1 Waste treatment methods	
13.1. Waste treatment methods	May be vented to atmosphere in a well ventilated place. Discharge to atmosphere in large quantities should be avoided. Do not discharge into any place where its accumulation could be dangerous. Return unused product in original cylinder to supplier.
13.2. Additional information	: External treatment and disposal of waste should comply with applicable local and/or national regulations
SECTION 14: Transport informatio	n
14.1. UN number	
UN-No.	: 1013
14.2. UN proper shipping name	
Transport by road/rail (ADG)	: CARBON DIOXIDE
Transport by air (ICAO-TI / IATA-DGR)	: Carbon dioxide
Transport by sea (IMDG)	: CARBON DIOXIDE



SDS Ref.: ALH062


14.3. Transport hazard class(es)	
Labelling	
	2
	2.2 : Non-flammable, non-toxic gases
Transport by road/rail (ADG)	
Class	: 2
Hazchemcode	: 2T
Hazard identification number	: 20
Tunnel Restriction	: C/E - Tank carriage : Passage forbidden through tunnels of category C, D and E. Other carriage : Passage forbidden through tunnels of category E
Transport by air (ICAO-TI / IATA-DGR)	
Class / Div. (Sub. risk(s))	: 2.2
Transport by sea (IMDG)	
Class / Div. (Sub. risk(s))	: 2.2
Emergency Schedule (EmS) - Fire	: F-C
Emergency Schedule (EmS) - Spillage	: S-V
14.4. Packing group	
Transport by road/rail (ADR/RID)	: Not applicable
Transport by air (ICAO-TI / IATA-DGR)	: Not applicable
Transport by sea (IMDG)	: Not applicable
14.5. Environmental hazards	
Transport by road/rail (ADR/RID)	: None.
Transport by air (ICAO-TI / IATA-DGR)	: None.
Transport by sea (IMDG)	: None.
14.6. Special precautions for user	
Packing Instruction(s)	
Transport by road/rail (ADR/RID)	: P200
Transport by air (ICAO-TI / IATA-DGR)	
Passenger and Cargo Aircraft	: 200

Cargo Aircraft only

Transport by sea (IMDG)

: 200

: P200

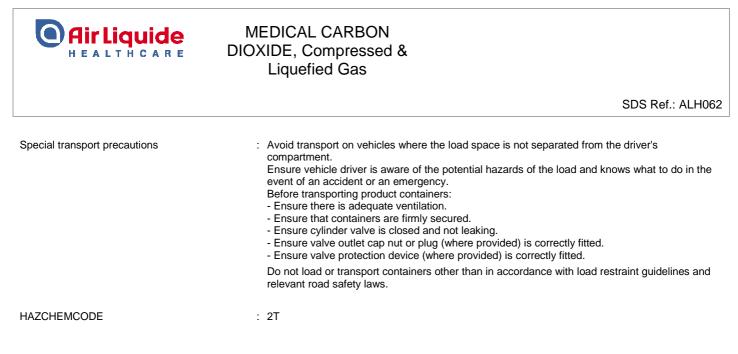

14.7. Transport in bulk according to Annex II of MARPOL 73/78 and the IBC Code

: Not applicable.

SECTION 15: Regulatory information		
15.1. Safety, health and environmental regulations/legislation specific for the substance or mixture		
National regulations		
Therapeutic Goods	: Registered presctiption medicine on the Australian Register of Therapeutic Goods (ARTG ID 32742)	
15.2. Chemical safety assessment		
	: A CSA does not need to be carried out for this product.	
SECTION 16: Other information		
Indication of changes	: Revised safety data sheet in accordance with commission regulation (EU) No 453/2010 and relevant Work, Health and Safety regulations.	
Abbreviations and acronyms	ATE - Acute Toxicity Estimate. CLP - Classification Labelling Packaging Regulation; Regulation (EC) No 1272/2008. REACH - Registration, Evaluation, Authorisation and Restriction of Chemicals Regulation (EC) No 1907/2006. EINECS - European Inventory of Existing Commercial Chemical Substances. CAS# - Chemical Abstract Service number. PPE - Personal Protection Equipment. LC50 - Lethal Concentration to 50 % of a test population. RMM - Risk Management Measures. PBT - Persistent, Bioaccumulative and Toxic. vPvB - Very Persistent and Very Bioaccumulative. STOT- SE : Specific Target Organ Toxicity - Single Exposure. CSA - Chemical Safety Assessment. EN - European Standard. UN - United Nations. ADR - European Agreement concerning the International Carriage of Dangerous Goods by Road. IATA - International Air Transport Association. IMDG code - International Maritime Dangerous Goods. RID - Regulations concerning the International Carriage of Dangerous Goods by Rail. WGK - Water Hazard Class. STOT - RE : Specific Target Organ Toxicity - Repeated Exposure.	
Training advice	: The hazard of asphyxiation is often overlooked and must be stressed during operator training.	

Full text of H-statements

Press. Gas (Liq.)	Gases under pressure : Liquefied gas
H280	Contains gas under pressure; may explode if heated.



SDS Ref.: ALH062

DISCLAIMER OF LIABILITY

: Before using this product in any new process or experiment, a thorough material compatibility and safety study should be carried out. Details given in this document are believed to be correct at the time of going to press. Whilst proper care has been taken in the preparation of this document, no liability for injury or damage resulting from its use can be accepted.